720

2022.10.15--388x128mm

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Operation Program of Lower Limbs

Prone position

Small flat head vertically downward

Hamstring	
Biceps femoris	20s-30sx3 times
Semitendinosus	20s-30sx3 times
Semimembranosus	20s-30sx3 times

The movement direction is from up to down.

Small flat head outward 45°

Gracilis	10s-15sx2 times
Adductor longus	10s-15sx2 times
Adductor brevis	10s-15sx2 times
Triceps surae	
Gastrocnemius	20s-30sx3 times
Soleus	20s-30sx3 times

The movement direction is from up to down.

Small round head vertically downward

Popliteus	
Knee straightening position	15s-30s
Knee bending position	15s-30s

Supine position

Small flat head vertically downward

Quadriceps femoris	
Rectus femoris	30s-45sx2 times
Vastus lateralis	30s-45sx2 times
Vastus medialis	30s-45sx2 times
Vastus intermedius	30s-45sx2 times

The movement direction is from up to down.

Small flat head vertically downward

Sartorius	30S-45S
Tensor fasciae latae	15s-20sx2 times
Tibialis anterior	15s-20sx2 times
Peroneus longus	15s-20sx2 times
Peroneus brevis	15s-20sx2 times

The movement direction is from up to down.

Use Method

Gonarthrits

Supine position

Small flat head vertically downward

Quadriceps femoris	
Rectus femoris	30s-45sx2 times
Vastus lateralis	30s-45sx2 times
Vastus medialis	30s-45sx2 times
Vastus intermedius	30s-45sx2 times


The movement direction is from up to down.

The patient slowly and slightly bends the knee.

Prone position

Small round head vertically downward

Popliteus	30s-45s
Knee straightening position	15s-30s
Knee bending position	15s-30s
Triceps surae	
Gastrocnemius	20s-30sx3 times
Soleus	20s-30sx3 times

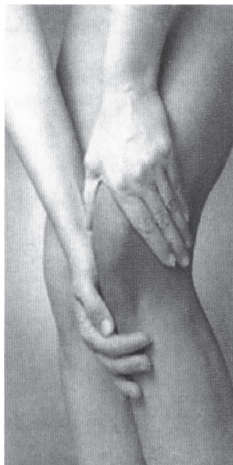Small round head vertically downward

Iliac bone edge - iliac bone anterior	30s-1min
Rectus femoris, inguen	1min
Tensor fasciae latae	15s-20sx2 times
Patellar ligament area	20sx3 times

Plantar fasciitis

Prone position (single side)

Small flat head vertically downward:

Hamstring	
Biceps femoris	20s-30sx3 times
Semitendinosus	20s-30sx3 times
Semimembranosus	20s-30sx3 times

The movement direction is from up to down.

U head vertically downward

Popliteus	
Knee straightening position	15s-30s
Knee bending position	15s-30s
Triceps surae	
Gastrocnemius	20s-30sx3 times
Soleus	20s-30sx3 times

The front foot sole touches the ground.

Small round head vertical to planta pedis:
Plantar fascia

Operation Program of Sacroiliac Part

Prone position (single side)

Small flat head vertically downward:

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

The direction is from inside up to outside down, along the muscular direction.

Hamstring	
Biceps femoris	20s-30sx3 times
Semitendinosus	20s-30sx3 times
Semimembranosus	20s-30sx3 times
From up to down	

Supine position

Tibialis anterior	15s-20sx2 times
Peroneus longus	15s-20sx2 times
Peroneus brevis	15s-20sx2 times

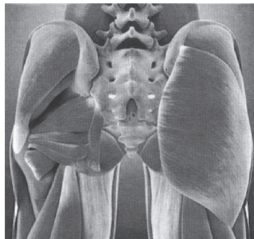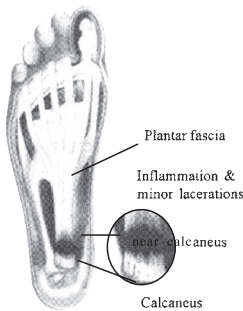The movement direction is from up to down.

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